

# SCREENi®

Viewing solution for Broncoflex®




- ▶ HIGH RESOLUTION SCREEN
- ▶ COMPACT AND MOBILE SYSTEM
- ▶ EASY TO USE AND INTUITIVE
- ▶ EXAMINATION MANAGEMENT INTERFACE
- ▶ EXPORT PHOTOS AND VIDEOS VIA USB

**axessvision**  
part of The Surgical Company

HIGH QUALITY VISION & SAFETY.  
ANYTIME, ANYWHERE.

# SCREENi® Viewing solution for Broncoflex®

High resolution touch panel

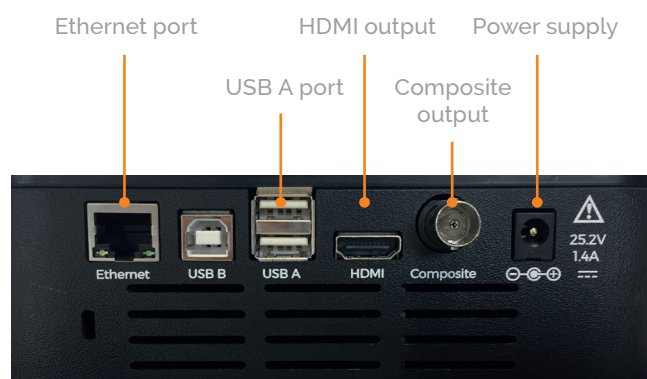

## HIGH RESOLUTION AND MOBILITY

**HIGH RESOLUTION SCREEN** 1280x800 with a bright touch panel offering optimal colour reproduction and contrast.

**COMPACT AND MOBILE SYSTEM** easily transportable over examination sites, the display may also be exported to an already existing medical monitor.

**SIMPLE TO USE AND INTUITIVE**, starts in seconds and immediately displays the video from the connected device. Videos and photos can be recorded during examination.

**EXAMINATION MANAGEMENT INTERFACE** making it possible to view photos and videos of performed exams.

**EXPORT RECORDED PHOTOS AND VIDEOS VIA USB** for inclusion in the patient file.

## ACCESSORIES



Robust mobile cart designed to connect to a SCREENi® tablet and fitted with provisions on which to hang Broncoflex® bronchoscopes in their sterile packaging.

Réf. 006.0003



Universal clamp for attaching SCREENi® to rails and tubes.

Réf. 006.0020

## TECHNICAL SPECIFICATIONS

ELECTRICAL	Power supply	External mains AC-DC Adapter 100 - 240V - 0,5 A 50 - 60 Hz
	Type of protection against electrocution	Class 2
	Battery type	Lithium-Ion battery
	Autonomy	3h30
	Charging time	2h
TOUCH SCREEN	Screen size	10 inches (21.6 x 13.5 cm)
	Resolution	1280 x 800
MISCELLANEOUS	IP Code	IP 30
	Dimensions (L x H x D)	310 x 220 x 70 mm
	Weight	2 Kg
	Accessories delivered with SCREENi®	Charger
	Processor	Embedded smart system
	Physical memory size	16 GB
	Input connector	Proprietary video connector 2 x USB Type A
COMPLIANCE	Output connector	HDMI
	MDD classification	Class IIa
ORDERING INFORMATION	SCREENi® 10" (including power cord and battery charger)	005.0008